# AWARENESS TECHNOLOGY, INC.

Founded in 1982, Awareness Technology has been supplying the international IVD market with high quality bench-top analyzers and reagents for biochemistry, ELISA, CLIA and ISE. Products include semi and fully automatic designs with famous brand names such as StatFax, ChroMate, LuMate, ChemWell, IONyte, and more.

Instruments are manufactured in the USA under both ISO and FDA Quality Systems.

With more than 100,000 instruments in the field, Awareness Technology is ready to supply customers with technical support including applications and service.

Our goal is reliable results at affordable prices.

# DELIVER RESULTS WITH CONFIDENCE

Designed and manufactured in the USA by

AWARENESS TECHNOLOGY, INC.

PO Box 1679, Palm City, FL 34991 USA

tel. (772) 283 6540 fax (772) 283 8020

web site: http://www.awaretech.com

email: info@awaretech.com

Available worldwide.

For sales information info@awaretech.com

For technical support cw-support@awaretech.com



Automatically process

ELISA and CLIA

in the same plate

in the same run

Patented dual-function reader automatically switches between absorbance and luminescence









#### **USER-PROGRAMMABLE OPEN SYSTEM**

Utilize a wide range of glow-luminescent assays and standard colorimetric ELISA tests including hormones, infectious diseases, tumor markers, auto-immune, and more.

# Perform automated ELISA and CLIA simultaneously

#### **IMPROVE QUALITY OF CARE**

- Easy to operate
- High-precision pipetting
- Software prevents timing conflicts
- Complete QC package

#### **COST EFFECTIVE OPERATION**

- Use reagents of your choice
- Consolidate two instruments in one
- No proprietary consumables
- Compact footprint saves bench space
- Robust design
- No disposable tips needed


#### **DATA MANAGEMENT**

- Password security
- Patient data base software included
- Unlimited number of customized protocols
- LIS integration
- Flexible reporting styles


#### **ROBUST AND RELIABLE**

- Self-monitoring of performance
- Simple economical maintenance
- Built in safety features
- Long-life IAD filters
- ISO13485, CE Marked

EASY
RELIABLE
ECONOMICAL

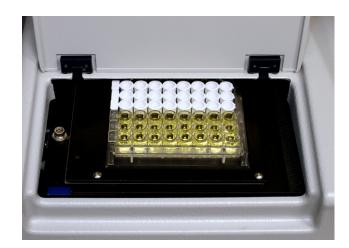



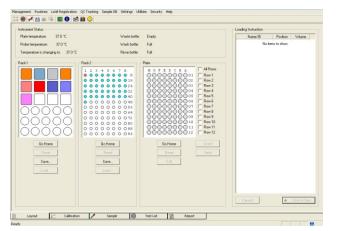

#### **SMART SOFTWARE**

The software automatically manages each assay separately so that multiple assays can be processed on the same plate even if they have different timings, reagents, and reading modes.





# **CAPACITY**

Samples: removable rack holds 96 serum vials or 12x75mm tubes, pre-dilution on-board the standard rack holds 27 various size bottles, optional alternate racks available standard 96-well microplates or strip trays - clear for ELISA and opaque white for

Chemiluminescent

## **OPTICAL PERFORMANCE**

**Absorbance** 

Linear Range: 0.0 to 3.0A Filters: long-life IAD

Wavelengths: 405, 450, 492, and 630nm

Luminescence

Sensitivity/Detection Limit: ALP 1x10<sup>-21</sup> moles HRP 1x10<sup>-18</sup> moles

Linear range: to 10<sup>6</sup> RLU
Cross talk: Less than 2.5 x 10<sup>-6</sup>

# **LIQUID HANDLING**

Volumes: from 2uL to 1.95 mL with precision better than 2%CV, single and multi dispense

Probe: detects liquid surface, automatic washing and rinsing Strip Washer: 8-probe manifold, two wash buffers on board

#### **CALCULATION MODES**

Single and multi-calibrator Linear and point-to-point Linear and non-regressions Log-logit, cubic spline Curve storage
Curve editing
Levy-Jennings plots
Cutoff calculations

### **PHYSICAL**

Dimensions: 86 cm W x 51 cm L x 40 cm H, approx. weight = 35 kg34 "W x 20 "L x 16 "H, approx. weight = 77 lbs.

Power: 100 - 250 VAC, 50-60 Hz

USB port and cable supplied

ChemWell® Fusion is supplied only by qualified, factory-trained distributors.